




CARFIL QUALITY
YOUR SUPPLIER FOR ALL LABORATORY NEEDS

Data Loggers Mercury Software

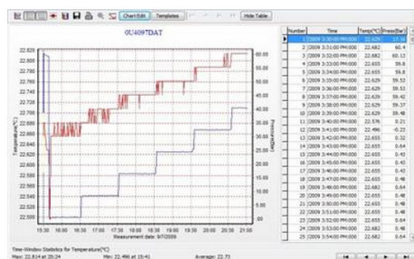

Overview

Mercury is our graphical Windows supporting software for the Pharmaceutical Industry. It is suitable for DST micro-T and DST nano-T. Just like our other software, Mercury offers various options for measurement setup and graphical and tabular view of retrieved data.

Description

In Mercury, the user sets the start time, start date and sampling interval before starting the recorder. A single interval can be defined, in second(s), minute(s) or hour(s).

Up to seven different intervals can be set for the same measurement sequence. Number of measurements is defined for each interval. This is especially useful when more frequent measurements are needed at certain time periods.

After recovering the DST, recorded data is uploaded to the software where the results are displayed in both graphic and tabular form.



PAVAN SERVICE BVBA

Beyntellus 3 - B-2360 Oud-Turnhout - Tel. +32 (0)14 45 13 10 - admin@carfil.be - BTW BE 0420.681.872

www.carfil.be